

**NIST Special Publication 260**  
**NIST SP 260-218r1**

# **Value Assignment of Reference Material 8404**

*Almond Flour for Allergen Detection*

Melissa M. Phillips  
David M. Bunk  
Ashley B. Green  
James H. Yen

This publication is available free of charge from:  
<https://doi.org/10.6028/NIST.SP.260-218r1>

**NIST Special Publication 260**  
**NIST SP 260-218r1**

# **Value Assignment of Reference Material 8404**

*Almond Flour for Allergen Detection*

Melissa M. Phillips  
*Chemical Sciences Division  
Material Measurement Laboratory*

David M. Bunk  
Ashley B. Green  
*Biomolecular Measurement Division  
Material Measurement Laboratory*

James H. Yen  
*Statistical Engineering Division  
Information Technology Laboratory*

This publication is available free of charge from:  
<https://doi.org/10.6028/NIST.SP.260-218r1>

October 2026

INCLUDES UPDATES AS OF 10-07-2026; SEE APPENDIX A



U.S. Department of Commerce  
*Howard Lutnick, Secretary*

National Institute of Standards and Technology  
*Arvind Raman, NIST Director and Under Secretary of Commerce for Standards and Technology*

Certain equipment, instruments, software, or materials, commercial or non-commercial, are identified in this paper in order to specify the experimental procedure adequately. Such identification does not imply recommendation or endorsement of any product or service by NIST, nor does it imply that the materials or equipment identified are necessarily the best available for the purpose.

#### **NIST Technical Series Policies**

[Copyright, Use, and Licensing Statements](#)

[NIST Technical Series Publication Identifier Syntax](#)

#### **Publication History**

Approved by the NIST Editorial Review Board on 2026-10-06

Supersedes NIST Special Publication 260-218 (December 2021) <https://doi.org/10.6028/NIST.SP.260-218>

#### **How to Cite this NIST Technical Series Publication**

Phillips MM, Bunk DM, Green AB, Yen JH (2026) Value Assignment of Reference Material 8404 Almond Flour for Allergen Detection. (National Institute of Standards and Technology, Gaithersburg, MD), NIST Special Publication (SP) 260-218r1. <https://doi.org/10.6028/NIST.SP.260-218r1>

#### **Author ORCID iDs**

Melissa M. Phillips: 0000-0003-0477-7637

David M. Bunk: 0009-0007-2187-4224

Ashley B. Green: 0000-0002-2065-4218

James H. Yen: 0000-0001-5242-7629

## **Abstract**

The National Institute of Standards and Technology (NIST) recently released RM 8404 Almond Flour for Allergen Detection, which is intended for harmonizing measurements of allergenic proteins in foods. The material was purchased from a commercial vendor and data was obtained from an interlaboratory comparison exercise and collaborating laboratories. A description of the material, results, and data analysis is discussed in the following report.

## **Keywords**

Allergen; Almond; Protein; Proximates; Reference Material.

## Table of Contents

|                                                                                                  |           |
|--------------------------------------------------------------------------------------------------|-----------|
| <b>1. Introduction.....</b>                                                                      | <b>1</b>  |
| <b>2. Material.....</b>                                                                          | <b>2</b>  |
| 2.1. Acquisition & Packaging.....                                                                | 2         |
| 2.2. Storage.....                                                                                | 5         |
| <b>3. Experimental Procedures.....</b>                                                           | <b>6</b>  |
| 3.1. Interlaboratory Studies for Value Assignment.....                                           | 6         |
| 3.2. Collaborating Laboratories for Value Assignment.....                                        | 6         |
| 3.3. Statistical Approaches for Value Assignment.....                                            | 7         |
| 3.4. Screening for Trace Allergen Contaminants.....                                              | 7         |
| 3.4.1. Eurofins GeneScan.....                                                                    | 7         |
| 3.4.2. Food Allergen Research and Resource Program (FARRP) at the University of Nebraska Lincoln | 8         |
| 3.4.3. Neogen Veratox Hazelnut ELISA.....                                                        | 8         |
| 3.4.4. Neogen Veratox Peanut ELISA.....                                                          | 8         |
| 3.4.5. BioFront Technologies MonoTrace for Pecan ELISA and Walnut ELISA.....                     | 9         |
| 3.4.6. 3M Macadamia Protein ELISA.....                                                           | 9         |
| 3.4.7. Microbac Laboratories.....                                                                | 10        |
| 3.5. Customer Feedback.....                                                                      | 10        |
| 3.5.1. Bia Diagnostics, LLC.....                                                                 | 10        |
| 3.5.2. Eurofins Analytik GmbH.....                                                               | 10        |
| 3.5.3. Eurofins Immunolab.....                                                                   | 11        |
| 3.5.4. Hygenia.....                                                                              | 11        |
| 3.5.5. Neogen.....                                                                               | 11        |
| <b>4. Results and Discussion.....</b>                                                            | <b>12</b> |
| 4.1. Proximates.....                                                                             | 12        |
| 4.1.1. Fat.....                                                                                  | 12        |
| 4.1.2. Protein.....                                                                              | 13        |
| 4.1.3. Carbohydrates.....                                                                        | 15        |
| 4.1.4. Ash.....                                                                                  | 15        |
| 4.1.5. Solids.....                                                                               | 16        |
| 4.1.6. Calories.....                                                                             | 16        |
| 4.2. Trace Allergen Contaminants.....                                                            | 17        |
| 4.2.1. Eurofins GeneScan.....                                                                    | 17        |
| 4.2.2. Food Allergen Research and Resource Program (FARRP) at the University of Nebraska Lincoln |           |

|                                           |           |
|-------------------------------------------|-----------|
| 4.2.3. Microbac Laboratories.....         | 18        |
| 4.3. Customer Feedback.....               | 18        |
| 4.3.1. Bia Diagnostics, LLC .....         | 18        |
| 4.3.2. Eurofins Analytik GmbH .....       | 19        |
| 4.3.3. Eurofins Immunolab .....           | 20        |
| 4.3.4. Hygenia.....                       | 20        |
| 4.3.5. Neogen .....                       | 20        |
| <b>5. Conclusions.....</b>                | <b>21</b> |
| 5.1. Value Assignment for Proximates..... | 21        |
| 5.2. Trace Allergen Contaminants.....     | 22        |
| 5.3. Customer Feedback.....               | 22        |
| <b>6. References .....</b>                | <b>23</b> |
| <b>Appendix A. Change Log.....</b>        | <b>25</b> |

## List of Tables

|                                                                                                                                                     |           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------|-----------|
| <b>Table 1. Results of August 2019 allergen screening by Romer Labs for candidate almond flour.....</b>                                             | <b>2</b>  |
| <b>Table 2. Results of September 2019 allergen screening by BioFront Laboratories for candidate almond flour. ....</b>                              | <b>3</b>  |
| <b>Table 3. Pre- and post-packaging homogeneity testing report from High-Purity Standards for RM 8404 Almond Flour for Allergen Detection. ....</b> | <b>5</b>  |
| <b>Table 4. Participants in the proximates study of HAMQAP Exercise 5.....</b>                                                                      | <b>6</b>  |
| <b>Table 5. Summary of Results for Fat, % .....</b>                                                                                                 | <b>12</b> |
| <b>Table 6. Summary of Eurofins Results for Protein, % .....</b>                                                                                    | <b>13</b> |
| <b>Table 7. Summary of HAMQAP Exercise 5 Results for Protein, %.....</b>                                                                            | <b>14</b> |
| <b>Table 8. Summary of Results for Carbohydrates, % .....</b>                                                                                       | <b>15</b> |
| <b>Table 9. Summary of Results for Ash, % .....</b>                                                                                                 | <b>15</b> |
| <b>Table 10. Summary of Results for Solids, % .....</b>                                                                                             | <b>16</b> |
| <b>Table 11. Summary of Results for Calories, kcal/100 g.....</b>                                                                                   | <b>16</b> |
| <b>Table 12. Results provided by Eurofins for hazelnut allergen in RM 8404 using R-Biopharm Hazelnut Allergen ELISA.....</b>                        | <b>17</b> |
| <b>Table 13. Results provided by FARRP for various allergens in RM 8404 .....</b>                                                                   | <b>17</b> |
| <b>Table 14. Results provided by Microbac Laboratories for various allergens in RM 8404 using R-Biopharm SureFood PCR Kits .....</b>                | <b>18</b> |
| <b>Table 15. Summary of cross reactivity and contamination results provided by Bia Diagnostics ELISA testing of RM 8404.....</b>                    | <b>18</b> |
| <b>Table 16. Results provided by Eurofins Analytik GmbH for RM 8404. ....</b>                                                                       | <b>19</b> |

|                                                                                                                           |           |
|---------------------------------------------------------------------------------------------------------------------------|-----------|
| <b>Table 17. Results for almond ELISA kit reactivity provided by Eurofins Immunolab for dilutions of RM 8404. ....</b>    | <b>20</b> |
| <b>Table 18. Results provided by Neogen for ELISA and lateral flow kit testing of materials spiked with RM 8404. ....</b> | <b>20</b> |
| <b>Table 19. Summary of Estimates for Proximates in RM 8404 .....</b>                                                     | <b>21</b> |
| <b>Table 20. Summary of all results for trace allergen contamination in RM 8404.....</b>                                  | <b>22</b> |

|                                                                                                                                                     |           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------|-----------|
| <b>Table 1. Results of August 2019 allergen screening by Romer Labs for candidate almond flour.....</b>                                             | <b>2</b>  |
| <b>Table 2. Results of September 2019 allergen screening by BioFront Laboratories for candidate almond flour. ....</b>                              | <b>3</b>  |
| <b>Table 3. Pre- and post-packaging homogeneity testing report from High-Purity Standards for RM 8404 Almond Flour for Allergen Detection. ....</b> | <b>5</b>  |
| <b>Table 4. Participants in the proximates study of HAMQAP Exercise 5.....</b>                                                                      | <b>6</b>  |
| <b>Table 5. Summary of Results for Fat, % .....</b>                                                                                                 | <b>12</b> |
| <b>Table 6. Summary of Eurofins Results for Protein, % .....</b>                                                                                    | <b>13</b> |
| <b>Table 7. Summary of HAMQAP Exercise 5 Results for Protein, %.....</b>                                                                            | <b>14</b> |
| <b>Table 8. Summary of Results for Carbohydrates, % .....</b>                                                                                       | <b>15</b> |
| <b>Table 9. Summary of Results for Ash, % .....</b>                                                                                                 | <b>15</b> |
| <b>Table 10. Summary of Results for Solids, % .....</b>                                                                                             | <b>16</b> |
| <b>Table 11. Summary of Results for Calories, kcal/100 g.....</b>                                                                                   | <b>16</b> |
| <b>Table 12. Results provided by Eurofins for hazelnut allergen in RM 8404 using R-Biopharm Hazelnut Allergen ELISA.....</b>                        | <b>17</b> |
| <b>Table 13. Results provided by FARRP for various allergens in RM 8404 .....</b>                                                                   | <b>17</b> |
| <b>Table 14. Results provided by Microbac Laboratories for various allergens in RM 8404 using R-Biopharm SureFood PCR Kits .....</b>                | <b>18</b> |
| <b>Table 15. Summary of cross reactivity and contamination results provided by Bia Diagnostics ELISA testing of RM 8404.....</b>                    | <b>18</b> |
| <b>Table 16. Results provided by Eurofins Analytik GmbH for RM 8404. ....</b>                                                                       | <b>19</b> |
| <b>Table 17. Results for almond ELISA kit reactivity provided by Eurofins Immunolab for dilutions of RM 8404. ....</b>                              | <b>20</b> |
| <b>Table 18. Results provided by Neogen for ELISA and lateral flow kit testing of materials spiked with RM 8404. ....</b>                           | <b>20</b> |
| <b>Table 19. Summary of Estimates for Proximates in RM 8404 .....</b>                                                                               | <b>21</b> |
| <b>Table 20. Summary of all results for trace allergen contamination in RM 8404.....</b>                                                            | <b>22</b> |

## List of Figures

|                                                                                                                                                       |           |
|-------------------------------------------------------------------------------------------------------------------------------------------------------|-----------|
| <b>Figure 1. Report of facility cleaning quality from High-Purity Standards prior to receipt of RM 8404 Almond Flour for Allergen Detection. ....</b> | <b>4</b>  |
| <b>Figure 2. Protein Mass Fraction Results Reported by Eurofins as a Function of Box Number .....</b>                                                 | <b>14</b> |

## **Introduction**

In 2017 and 2019, NIST held workshops to identify needs of the food industry and federal regulators. Among other things, NIST was asked to continue production of food-matrix SRMs for use by laboratories making measurements in support of food safety [1]. One described need was for additional commodity materials for method comparison and harmonization of food allergen testing, as laboratories need a means for demonstrating method validity and accuracy when analyzing food products. RM 8404 Almond Flour for Allergen Detection was also requested by various stakeholders through AOAC INTERNATIONAL to assist in the evaluation of allergen determination in food matrices. NIST currently offers commodity reference materials for many important allergenic foods (e.g., milk, egg, wheat, soy, peanuts, fish, shellfish), but RMs for tree nut allergens are currently unavailable from NIST or any other reference material producer. Availability of reference materials for tree nuts, such as almond, will facilitate the development and harmonization of methods for detecting trace levels of these allergenic foods in finished products.

## 1. Material

### 1.1. Acquisition & Packaging

Based on the intended use of this RM, selection of a material that contains proteins from a single nut source was critical. Numerous manufacturers were evaluated via websites and product claims and several contacted to determine whether the nut products that they produce are likely to be pure or may have come in contact with other tree nuts or allergenic foods. In August 2019, 2.2 kg (1 pound) of blanched almond flour was purchased from Mandelin (San Luis Obispo, CA). The candidate material was aliquoted into 100-g packets and sent in blind triplicate to Romer Labs (Union, MO) to be tested for the presence of protein allergens from cashew, wheat, hazelnut, macadamia, peanut, pecan, pine nut, pistachio, and walnut. Results of the allergen testing conducted by Romer Labs in August 2019 are described in Table 1, which identified possible presence of macadamia protein in the sample. After consulting with experts, the presence of macadamia protein in the almond flour was deemed questionable given the low prevalence of macadamia nut production in California and the potential for cross-reactivity of the enzyme-linked immunosorbent assay (ELISA) utilized by Romer Labs.

**Table 1. Results of August 2019 allergen screening by Romer Labs for candidate almond flour.**

| Parameter | Result<br>mg/kg (ppm) |      |      | Method                               | LOD <sup>a</sup> | LOQ <sup>b</sup> |
|-----------|-----------------------|------|------|--------------------------------------|------------------|------------------|
|           | A                     | B    | C    |                                      | mg/kg<br>(ppm)   | mg/kg<br>(ppm)   |
| Cashew    | <LOD                  | <LOD | <LOD | US-AgraQuant ELISA Allergen Test Kit | 0.2              | 2.0              |
| Gluten    | <LOD                  | <LOD | <LOD | US-AgraQuant ELISA Allergen Test Kit | 2.0              | 4.0              |
| Gluten    | <LOD                  | <LOD | <LOD | US-RIDASCREEN Gliadin Test Kit       | 3.0              | 5.0              |
| Hazelnut  | <LOD                  | <LOD | <LOD | US-AgraQuant ELISA Allergen Test Kit | 0.3              | 1.0              |
| Macadamia | 11.7                  | 12.7 | 10.5 | US-AgraQuant ELISA Allergen Test Kit | 1.0              | 2.0              |
| Peanut    | <LOD                  | <LOD | <LOD | US-AgraQuant ELISA Allergen Test Kit | 0.1              | 1.0              |
| Pecan     | <LOD                  | <LOD | <LOD | US-Pecan ELISA                       | 1.0              | 2.0              |
| Pine Nut  | <LOD                  | <LOD | <LOD | US-Pine Nut ELISA                    | 0.7              | 1.5              |
| Pistachio | <LOD                  | <LOD | <LOD | US-AgraQuant ELISA Allergen Test Kit | 0.1              | 1.0              |
| Walnut    | <LOD                  | <LOD | <LOD | US-AgraQuant ELISA Allergen Test Kit | 0.3              | 2.0              |

<sup>a</sup> LOD = limit of detection

<sup>b</sup> LOQ = limit of quantification

In September 2019, a 30-g aliquot of the candidate almond material was sent to BioFront Technologies (Tallahassee, FL) to be evaluated using different ELISAs than those used at Romer Labs. The material was tested for the presence of protein allergens from almond, brazil nut, cashew, hazelnut, macadamia nut, peanut, pecan, pine nut, pistachio, and walnut. Results of the allergen testing conducted by BioFront Technologies in September 2019 are described in Table 2. The candidate material was found to contain less than 1 mg/kg (ppm) macadamia protein, indicating that the positive result found by Romer Labs may have been a result of assay cross-reactivity.

**Table 2. Results of September 2019 allergen screening by BioFront Laboratories for candidate almond flour.**

| Parameter  | Result<br>mg/kg (ppm) | MonoTrace ELISA                 |                                 |                                  |
|------------|-----------------------|---------------------------------|---------------------------------|----------------------------------|
|            |                       | ROQ <sup>a</sup><br>mg/kg (ppm) | LOD <sup>b</sup><br>mg/kg (ppm) | LLOQ <sup>c</sup><br>mg/kg (ppm) |
| Almond     | >ULOQ <sup>d</sup>    | 1-40                            | 0.23                            | 1                                |
| Brazil Nut | <LLOQ                 | 1-40                            | 0.14                            | 1                                |
| Cashew     | <LLOQ                 | 1-40                            | 0.12                            | 1                                |
| Hazelnut   | <LLOQ                 | 1-40                            | 0.04                            | 1                                |
| Macadamia  | <LLOQ                 | 2-80                            | 0.13                            | 2                                |
| Peanut     | <LLOQ                 | 1-40                            | 0.24                            | 1                                |
| Pecan      | <LLOQ                 | 1-40                            | 0.17                            | 1                                |
| Pine Nut   | <LLOQ                 | 1-40                            | 0.24                            | 1                                |
| Pistachio  | <LLOQ                 | 1-40                            | 0.12                            | 1                                |
| Walnut     | <LLOQ                 | 1-40                            | 0.22                            | 1                                |

<sup>a</sup> ROQ = range of quantification

<sup>b</sup> LOD = limit of detection

<sup>c</sup> LLOQ = lower limit of quantification

<sup>d</sup> ULOQ = upper limit of quantification.

Final packaging of the RM was completed by High-Purity Standards (North Charleston, SC). Prior to receiving materials, Neogen Reveal kits for Multi-Treenut were shipped by NIST to High-Purity Standards. High-Purity Standards was asked to use these kits to test for cross-contact with other allergen-containing materials during packaging. In December 2019, the room and equipment at High-Purity Standards were thoroughly cleaned and the Reveal kit used to test swabs of the tabletop and packaging components. No tree nuts were detected prior to receipt of the almond flour, as shown in the report from High-Purity Standards (Figure 1).




December 30, 2019

National Institute of Standards and Technology  
100 Bureau Drive  
Bldg 301, Room B130  
Gaithersburg, MD 20899-1410

Reference: Allergen test report for contract 1333ND19PNB640946P20001 sales Order 17598

Dear Dr. Melissa Phillips,

High-Purity Standards received a test kit from NIST for the Neogen Reveal® for Multi-Treenut test. It was stored in a refrigerator maintained at 2°C – 8°C. Prior to the test, the room and equipment were thoroughly cleaned. The test was then performed according to manufacturer's instructions. We tested the tabletop and inner surfaces of the polybag, drum, air vent, balance, weigh-boat and gas nozzle. Tree nut was not detected on all equipment confirming that they are all clean and free of tree nut. Shown below is a representative image of the results for the tabletop.

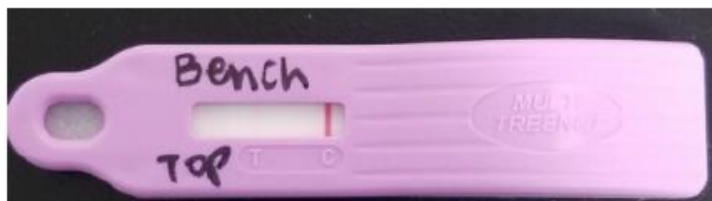

HPS is now ready to receive the almond powder.

**Participants:**

**Cleaning:** Carley Freeland – Inorganic Lab Technician  
**Analysis:** Tommy Breen – Inorganic Manufacturing Manager

Sincerely,

Tommy Breen, Inorganic Lab Manager



7221 Investment Drive • North Charleston, SC 29418  
843.767.7900 • info@highpuritystandards.com • www.highpuritystandards.com

Figure 1. Report of facility cleaning quality from High-Purity Standards prior to receipt of RM 8404 Almond Flour for Allergen Detection.

In January 2020, 34 kg (75 pounds) of blanched almond flour was purchased from Mandelin and shipped directly to High-Purity Standards. Prior to packaging, the material was blended for 8 h. To verify the homogeneity of the blended flour, six samples were analyzed by High-Purity Standards using inductively coupled plasma optical emission spectrometry (ICP-OES) following microwave digestion. The results for Ca, K, Mg, Mn, S, and Zn (Table 3) indicate that sufficient blending was achieved, and that the material could be packaged as requested. Aliquots ( $5.0 \text{ g} \pm 0.1 \text{ g}$ ) of the almond flour were weighed using a static free container on a platform balance with the accuracy of  $\pm 0.0001 \text{ g}$  and immediately transferred to 0.10 mm (4 mil or 0.004 in) plastic polyethylene bags through a solid funnel. The polyethylene bags were flushed with dry nitrogen and immediately heat sealed and over-packed aluminized plastic packets with two 0.5 g packets of SORB-IT while being purged with dry nitrogen before double sealing. The aluminized packets were placed in rows in 35.56 cm x 35.56 cm x 35.56 cm (14 in x 14 in x 14 in) cardboard boxes. The packets were arranged into sections (1 through 4), placed from back to front of the box. The front and back of each box were marked, and the boxes were numbered sequentially and sealed. A total of 10 boxes were produced containing 500 packets each, giving a grand total of 5000 packets. Fourteen packets were removed after packaging for additional homogeneity testing (Table 3). The remaining material was repackaged and included with the packaged material on a pallet and sealed in a plastic film for shipment to NIST. 10 boxes of packaged RM 8404 Almond Flour for Allergen Detection, as well as remaining unpackaged material, were received at NIST on February 25, 2020.

**Table 3. Pre- and post-packaging homogeneity testing report from High-Purity Standards for RM 8404 Almond Flour for Allergen Detection.**

| Analyte | % RSD                      |                             | Criteria (%) | Result |
|---------|----------------------------|-----------------------------|--------------|--------|
|         | Pre-Packaging <sup>a</sup> | Post-Packaging <sup>b</sup> |              |        |
| Ca      | 1.96                       | 2.23                        | $\leq 3.0$   | Pass   |
| K       | 1.94                       | 2.89                        | $\leq 3.0$   | Pass   |
| Mg      | 1.54                       | 1.91                        | $\leq 3.0$   | Pass   |
| Mn      | 1.91                       | 2.80                        | $\leq 3.0$   | Pass   |
| S       | 2.64                       | 2.58                        | $\leq 3.0$   | Pass   |
| Zn      | 1.69                       | 2.75                        | $\leq 3.0$   | Pass   |

<sup>a</sup> RSD for 6 randomly selected samples.

<sup>b</sup> RSD for 14 randomly selected samples.

## 1.2. Storage

The packets of RM 8404 have been stored at  $-20^\circ \text{C}$  at NIST since their receipt.

## 2. Experimental Procedures

### 2.1. Interlaboratory Studies for Value Assignment

RM 8404 was distributed in Exercise 5 of the NIST Health Assessment Measurements Quality Assurance Program (HAMQAP). Laboratories participating in the HAMQAP Exercise were provided with 3 packets of RM 8404 and were asked to use their in-house analytical methods to determine the mass fraction (percent) of proximates (fat, protein, carbohydrates, solids, and ash) as well as calories (kcal/100 g) in each packet. The quantitative results from this study are reported here in full, and the full report from Exercise 5 is published elsewhere [2]. Results were reported by the participants listed in Table 4, using the methods described in Section 3.1. The reported results from each participating organization have been assigned an arbitrary numeric code.

**Table 4. Participants in the proximates study of HAMQAP Exercise 5.**

| <b>Company</b>                                | <b>Location</b> | <b>Country</b> |
|-----------------------------------------------|-----------------|----------------|
| Advanced Botanical Consulting & Testing, Inc. | Tustin, CA      | USA            |
| Analytical Resource Labs                      | Lehi, UT        | USA            |
| BI Nutraceuticals                             | Sparks, NV      | USA            |
| Exact Scientific Services, Inc.               | Ferndale, WA    | USA            |
| Intertek Champaign Laboratories               | Champaign, IL   | USA            |
| SORA Labs                                     | Forsyth, MO     | USA            |

### 2.2. Collaborating Laboratories for Value Assignment

Eurofins Food Chemistry Testing US (Madison, WI) was provided with 14 samples of RM 8404 for determination of total protein. To establish the repeatability of the laboratory's method, 5 samples were provided from the beginning of the production lot. Nine additional samples were provided from across the production lot to evaluate material homogeneity.

Total protein was determined in 0.2 g to 0.3 g samples of RM 8404 using Dumas combustion based on AOAC Official Method 968.06 Protein (Crude) in Animal Feed [3]. In summary, the samples were combusted at  $\geq 850$  °C and the nitrogen generated was carried by CO<sub>2</sub> for quantitation by thermal conductivity. The nitrogen content determined in the samples was converted to crude protein using a generic conversion factor of 6.25, common for legumes, corn, and many animal proteins [4]. Jones factors for nuts and seeds, however, are lower and should be used in the case of almond testing. The data provided by the collaborating laboratory has been adjusted using a more appropriate Jones factor for almond measurement (5.18).

$$\text{crude protein (\%)} = \text{nitrogen (\%)} \times \text{conversion factor}$$

## **2.3. Statistical Approaches for Value Assignment**

Statistical analysis was provided by the NIST Statistical Engineering Division (SED). Where more than one method is available for a measured analyte, the estimated value is a weighted mean of the method estimates available for this analyte. The weighted mean used is the Dersimonian-Laird estimate [5], the uncertainty of which is estimated using a bootstrap procedure based on a Gaussian random effects model for the between-method effects [6-9]. The uncertainty for protein also includes a relative uncertainty component of 2.335% to account for the reproducibility of the Dumas method used for its value assignment. If only one method is available for an analyte, then that method estimate is the analyte estimate.

Very marked differences are often observed between the results from the different laboratories participating in an interlaboratory study. For the interlaboratory study, the estimate for each analyte is the weighted median of the individual laboratory means for that analyte, where the weights are based on a Laplace random effects model [10]. For this RM, the weighted median is equal to the unweighted median of laboratory means for all analytes. The uncertainty of the weighted median is estimated using a bootstrap procedure based on a Laplace random effects model for the between-laboratory and within-laboratory effects [6-10].

Some measurements from the interlaboratory studies were flagged by the analysts and omitted from the calculations. The deviance of these measurements from the others exceeded the usual variation, often differing by an order of magnitude or more. Other measurements may be questionable but could not be determined to be unrepresentative extreme outliers because of the sparseness and variation of the rest of the data. Some measurements were corrected for incorrect reporting units or incorrect Jones factors (for protein) and are noted in the sections below.

Some of the estimates and uncertainties in this report are purposely listed with more significant digits than is scientifically warranted. The relevant technical experts trim any estimates and uncertainties to the number of significant digits that are scientifically warranted prior to inclusion on the Reference Material Information Sheet as non-certified values [11].

## **2.4. Screening for Trace Allergen Contaminants**

### **2.4.1. Eurofins GeneScan**

Eurofins GeneScan (New Orleans, LA) was provided with 3 samples of RM 8404 for testing by R-Biopharm R6802 Hazelnut Allergen (ELISA). The RIDASCREEN FAST Hazelnut (Product R6802) is a sandwich enzyme immunoassay for the quantitative analysis of hazelnut (or hazelnut proteins) in food with a limit of quantification (LOQ) of 2.5 mg/kg (ppm) hazelnut [12]. The test principle is described below in an excerpt from reference [12].

The basis of the test is the antigen-antibody reaction. The wells of the microtiter strips are coated with specific antibodies against hazelnut proteins. By adding the standard or sample solution to the wells, present hazelnut protein will bind to the specific capture antibodies. The result is an antibody-antigen-complex. Components not bound by the

antibodies are then removed in a washing step. Then, antibody conjugated to peroxidase is added. This conjugate is bound to the Ab-Ag-complex. An antibody-antigen-antibody (sandwich) complex is formed. Any unbound conjugate is then removed in a washing step. Substrate/Chromogen are added to the wells and incubated. Bound conjugate converts the colorless chromogen into a blue product. The addition of the stop solution leads to a color change from blue to yellow. The measurement is made photometrically at 450 nm. The absorbance is proportional to the hazelnut content of the sample. The result is expressed in mg/kg hazelnut.

#### **2.4.2. Food Allergen Research and Resource Program (FARRP) at the University of Nebraska Lincoln**

The Food Allergen Research and Resource Program (FARRP) at the University of Nebraska Lincoln (Lincoln, NE) was provided with samples of RM 8404 for testing by numerous commercial approaches as described in the sections below. FARRP was provided with 3 samples of RM 8404 for each test.

#### **2.4.3. Neogen Veratox Hazelnut ELISA**

The Neogen Veratox for Hazelnut Allergen Quantitative Test (Product 8420) is a sandwich enzyme-linked immunosorbent assay intended for the full quantitative analysis or simple screening of hazelnut protein residues in food products with an LOQ of 2.5 mg/kg (ppm) total hazelnut [13]. The test principle is described below in an excerpt from reference [13].

The Veratox for Hazelnut Allergen test is a sandwich enzyme-linked immunosorbent assay (S-ELISA). Hazelnut protein is extracted from samples with a phosphate buffered salt solution (PBS) by shaking in a heated water bath, followed by centrifugation or filtration. Extracted hazelnut protein is sampled and added to antibody-coated wells (capture antibody) where it binds to the antibody during an incubation. Any unbound hazelnut protein is washed away and a second antibody (detector antibody), which is enzyme labeled, is added. The detector antibody binds to the already bound hazelnut protein. After a second wash, substrate is added. Color develops as a result of the presence of bound detector antibody. Red Stop reagent is added and the color of the resulting solution is observed. The test is read in a microwell reader to yield optical densities. The optical densities of the controls form a standard curve, and the sample optical densities are plotted against the curve to calculate the exact concentration of hazelnut expressed as parts per million of hazelnut.

#### **2.4.4. Neogen Veratox Peanut ELISA**

The Neogen Veratox for Peanut Allergen Quantitative Test (Product 8420) is a sandwich enzyme-linked immunosorbent assay intended for the quantitative analysis of peanut protein in food products with an LOQ of 2.5 mg/kg (ppm) total peanut [14]. The test principle is described below in an excerpt from reference [14].

The Veratox for Peanut Allergen test is a sandwich enzyme-linked immunosorbent assay (S-ELISA). Peanut protein residue is extracted from samples with a phosphate buffered salt solution (PBS) by shaking in a heated water bath, followed by centrifugation or filtration. Extracted peanut residue is sampled and added to antibody-coated wells (capture antibody) where it binds to the antibody during an incubation. Any unbound peanut residue is washed away and a second antibody (detector antibody), which is enzyme labeled, is added. The detector antibody binds to the already bound peanut residue. After a second wash, substrate is added. Color develops as a result of the presence of bound detector antibody. Red Stop reagent is added and the color of the resulting solution is observed. The test is read in a microwell reader to yield optical densities. The optical densities of the controls form a standard curve, and the sample optical densities are plotted against the curve to calculate the exact concentration of peanut protein.

#### **2.4.5. BioFront Technologies MonoTrace for Pecan ELISA and Walnut ELISA**

The BioFront Technologies MonoTrace for Pecan ELISA (Product PC4-EK-48) is monoclonal antibody-based assay intended for the qualitative or quantitative detection of pecan protein in food products with an LOQ of 1.0 mg/kg (ppm) total pecan [15]. BioFront Technologies MonoTrace for Walnut ELISA (Product WJ4-EK-48) is monoclonal antibody-based assay intended for the qualitative or quantitative detection of walnut protein in food products with an LOQ of 1.0 mg/kg (ppm) total walnut [15]. As summarized from the graphical representation in [15], the target protein residue is extracted from samples with a buffer solution. Extracted residue is sampled and added to monoclonal antibody-coated wells (capture antibody) where it binds to the antibody during an incubation. Any unbound residue is washed away in a series of wash steps and a second horseradish peroxidase (HRP)-conjugated monoclonal antibody (detector antibody) is added. The detector antibody binds to the already bound residue. After a second series of washes, substrate is added. Color develops as a result of the presence of bound detector antibody. HRP-quench solution is added to stop the reaction and the color of the resulting solution is observed. The test is read in a microwell reader to yield optical densities. The optical densities of the controls form a standard curve, and the sample optical densities are plotted against the curve to calculate the exact concentration of trace protein.

#### **2.4.6. 3M Macadamia Protein ELISA**

The 3M Macadamia Protein ELISA (Product E96MAC) is a sandwich enzyme-linked immunosorbent assay intended for the detection of macadamia proteins in clean-in-place water (CIP) final rinse water, environmental swab samples, food ingredients, and processed food products with an LOQ of 0.33 mg/kg (ppm) macadamia protein [16]. The test principle is described below in an excerpt from reference [16].

The 3M Macadamia Protein ELISA Kit utilizes a sandwich ELISA. The macadamia proteins present in the sample react with the anti-macadamia antibody, which have been adsorbed to the surface of polystyrene microtiter wells. After the removal of unbound

proteins by washing, anti-macadamia antibodies conjugated with horseradish peroxidase (HRP) are added. These enzyme-labeled antibodies form complexes with the previously bound macadamia protein. Following a second washing step, the enzyme bound to the immunosorbent is detected by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The color development from this enzymatic reaction varies directly with the concentration of macadamia protein in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of macadamia protein in the test sample. The quantity of macadamia protein in the test sample can be extrapolated from the standard curve, constructed from standards of known concentration, and adjusted to consider the sample dilution.

#### **2.4.7. Microbac Laboratories**

Microbac Laboratories (Oak Ridge, TN) was provided with 6 sets each of 3 samples of RM 8404 for testing by R-Biopharm SureFood Allergen PCR for peanut (S3603), hazelnut (S3602), cashew (S3615), pecan (S3618), pistachio (S3614), and walnut (S3607). SureFood Allergen PCR is a polymerase chain reaction technology that utilizes an internal amplification control for qualitative and quantitative detection of allergenic foods with a limit of detection of 0.4 mg/kg (ppm) of DNA for peanut and all tree nuts of interest [17]. As summarized in reference [17], the substance to be tested is lysed in buffer and proteinase K at 65 °C for one hour. After centrifugation and filtration via a spin filter, the DNA is bound to a spin filter, washed several times with wash buffer and eluted with of elution buffer. In a two-step thermal profile, the DNA is amplified for 45 cycles. A positive result for a qualitative test shows an exponential curve and a cycle threshold (Ct) value.

### **2.5. Customer Feedback**

Samples of RM 8404 Almond Flour for Allergen Detection (2 packets of 5 g of material, for a total of 10 g of material) were provided to interested stakeholders for evaluation of fitness-for-purpose. Several laboratories provided feedback to NIST regarding their use of RM 8404, as summarized below.

#### **2.5.1. Bia Diagnostics, LLC**

Bia Diagnostics, LLC (Colchester, VT) utilized 3M Protein ELISA Kits for coconut, almond, hazelnut, pecan, macadamia, walnut, pistachio, cashew, and Brazil nut testing, which are quantitative kits utilizing polyclonal antibodies.

#### **2.5.2. Eurofins Analytik GmbH**

Eurofins Analytik GmbH (Hamburg, Germany) utilized internal methods for allergen testing, including ELISA testing for  $\beta$ -lactoglobulin and casein (milk proteins), cashew, peanut, coconut, hazelnut, lupin, macadamia, Brazil nut, pecan, sesame, soy, egg, walnut, and mustard, and polymerase chain reaction (PCR) testing for cashew, peanut, fish, oat, hazelnut, lupin, pistachio,

celery, mustard, sesame, soy, walnut, and wheat. All tests were conducted in duplicate, and the results of the single determinations did not vary significantly. Additionally, paprika powder and chocolate were spiked with RM 8404 leading to a final concentration of 1 % (10000 mg/kg, 10000 ppm) almond and tested with the almond ELISA.

### **2.5.3. Eurofins Immunolab**

Eurofins Immunolab (Kassel, Germany) utilized internal ELISA methods (three kit lots) for allergen testing. Proteins from RM 8404 were extracted according to the kit instructions for use and diluted in extraction buffer.

### **2.5.4. Hygenia**

Hygenia (Camas, Sevilla, Spain) used RM 8404 in quality control and cross reactivity studies. No additional information was provided.

### **2.5.5. Neogen**

Neogen (Lansing, MI) used Veratox Almond (8440) kits to evaluate the samples via ELISA. 1000 µg/ml samples were prepared by adding 1 mg of the reference material to 1 mL of phosphate-buffered saline (PBS). The sample was vortexed until dissolved and then diluted 1:10 to prepare a 100 µg/mL spike stock. Five grams of material (rice or PBS) was spiked at (0, 5, 10, and 20) mg/kg (ppm) commodity using the 100 µg/mL stock. The samples were extracted and tested in duplicate following the kit insert.

Neogen also used their lateral flow assay Reveal 3D Almond (902086G) to evaluate RM 8404. The 100 µg/mL spike stock described above was used to spike 0.25 g of commodity (rice flour and PBS) at (0, 5, and 10) µg/mL. These samples, in addition to the 1000 µg/mL spike stock described above, were extracted and tested in duplicate via kit insert.

### 3. Results and Discussion

#### 3.1. Proximates

Results for proximates provided by Exercise 5 of HAMQAP and by Eurofins Food Chemistry Testing are summarized in the sections below. All results were provided on an as-received basis.

##### 3.1.1. Fat

The fat values reported by laboratories participating in HAMQAP Exercise 5 are summarized in Table 5. The table also provides several summary values: *N* = number of values, Mean = mean of values, SD = standard deviation of values, and Pooled SD = square-root of the sum of the squared within-participant standard deviations.

**Table 5. Summary of Results for Fat, %**

| Lab              | A     | B     | C     | Mean | SD  | Method                              |
|------------------|-------|-------|-------|------|-----|-------------------------------------|
| E002             | 58.64 | 58.26 | 59.01 | 58.6 | 0.4 | Sum of fatty acids as triglycerides |
| E030             | 51.5  | 50    | 49.7  | 50.4 | 1.0 | Sum of fatty acids as triglycerides |
| E033             | 62.34 | 59.6  | 59.2  | 60.4 | 1.7 | Sum of fatty acids as triglycerides |
| E047             | 53.71 | 53.09 |       | 53.4 | 0.4 | Rose-Gottlieb/Monjonier             |
| <i>N</i> :       |       |       |       | 4    |     |                                     |
| Mean, Pooled SD: |       |       |       | 55.9 | 2.0 |                                     |
| SD:              |       |       |       | 4.4  |     |                                     |


### 3.1.2. Protein

The protein values reported by Eurofins Food Chemistry Testing are summarized in Table 6. In the first column of data, as reported by Eurofins, the results for total nitrogen were converted to protein using a Jones factor of 6.25. In the second column of data, the results have been adjusted by NIST to reflect a proper Jones factor of 5.18. The table also provides several summary values: *N* = number of values, Mean = mean of values, and SD = standard deviation of values.

**Table 6. Summary of Eurofins Results for Protein, %**

| Box          | Jones Factor |      |
|--------------|--------------|------|
|              | 6.25         | 5.18 |
| 1-1          | 28.2         | 23.4 |
| 1-1A         | 28.1         | 23.3 |
| 1-1B         | 27.8         | 23.0 |
| 1-1C         | 27.8         | 23.0 |
| 1-1D         | 28.1         | 23.3 |
| 1-3          | 28.3         | 23.5 |
| 2-4          | 28.0         | 23.2 |
| 3-4          | 28.2         | 23.4 |
| 4-3          | 28.5         | 23.6 |
| 6-1          | 28.4         | 23.5 |
| 7-1          | 28.3         | 23.5 |
| 8-1          | 27.9         | 23.1 |
| 8-2          | 28.1         | 23.3 |
| 10-4         | 27.6         | 22.9 |
| <b>N:</b>    | 10           | 10   |
| <b>Mean:</b> | 28.1         | 23.3 |
| <b>SD:</b>   | 0.3          | 0.2  |

Figure 2 displays the protein results reported by Eurofins Food Chemistry Testing as a function of the sample box number using the Jones factor of 5.18. The blue circles in the figure represent the individual test results for each sample. The mean mass fraction is indicated by the solid line, and the dashed lines bound the interval  $\text{Mean} \pm 2 \times \text{SD}$ .

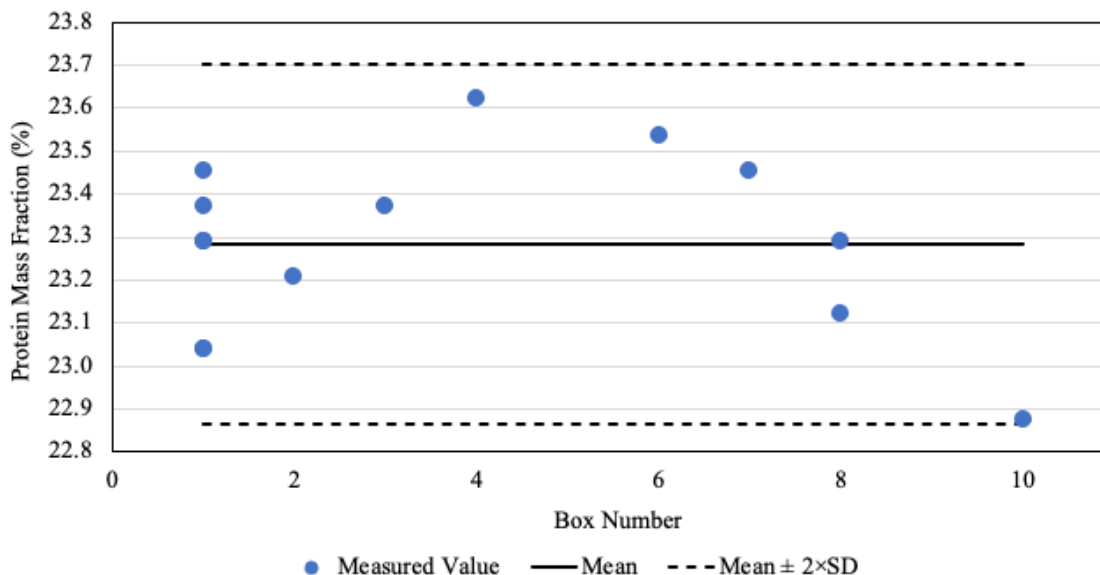

**Figure 2. Protein Mass Fraction Results Reported by Eurofins as a Function of Box Number**

The protein values reported by laboratories participating in HAMQAP Exercise 5 are summarized in Table 7. The upper portion of the table summarizes the results as reported by the participants in the study. The lower portion of the table repeats this information, replacing the results for two laboratories who reported using a Jones factor of 6.25. The table also provides several summary values:  $N$  = number of values, Mean = mean of values, SD = standard deviation of values, and Pooled SD = square-root of the sum of the squared within-participant standard deviations.

**Table 7. Summary of HAMQAP Exercise 5 Results for Protein, %**

| HAMQAP Exercise 5 (as reported by participants) |      |      |      |      |     |                       |
|-------------------------------------------------|------|------|------|------|-----|-----------------------|
| Lab                                             | A    | B    | C    | Mean | SD  | Method (Jones Factor) |
| E002                                            | 24.9 | 25.3 | 24.8 | 25   | 0.3 | Kjeldahl (5.18)       |
| E030                                            | 22.4 | 22.7 | 22.7 | 22.6 | 0.2 | Kjeldahl (5.18)       |
| E033                                            | 28.8 | 28.3 | 28.8 | 28.6 | 0.3 | Combustion (6.25)     |
| E047                                            | 28.0 | 28.1 |      | 28.1 | 0.1 | Combustion (6.25)     |
| <b>N:</b>                                       |      |      |      | 4    |     |                       |
| <b>Mean, Pooled SD:</b>                         |      |      |      | 25.9 | 0.4 |                       |
| <b>SD:</b>                                      |      |      |      | 2.6  |     |                       |

| HAMQAP Exercise 5 (corrected for Jones factor) |      |      |      |      |     |                       |
|------------------------------------------------|------|------|------|------|-----|-----------------------|
| Lab                                            | A    | B    | C    | Mean | SD  | Method (Jones Factor) |
| E002                                           | 24.9 | 25.3 | 24.8 | 25   | 0.3 | Kjeldahl (5.18)       |
| E030                                           | 22.4 | 22.7 | 22.7 | 22.6 | 0.2 | Kjeldahl (5.18)       |
| E033                                           | 23.9 | 23.5 | 23.9 | 23.7 | 0.2 | Combustion (5.18)     |
| E047                                           | 23.2 | 23.3 |      | 23.2 | 0.1 | Combustion (5.18)     |
| <b>N:</b>                                      |      |      |      | 4    |     |                       |
| <b>Mean, Pooled SD:</b>                        |      |      |      | 23.7 | 0.4 |                       |
| <b>SD:</b>                                     |      |      |      | 1.0  |     |                       |

### 3.1.3. Carbohydrates

The carbohydrate values reported by laboratories participating in HAMQAP Exercise 5 are summarized in Table 8. The table also provides several summary values: *N* = number of values, Mean = mean of values, SD = standard deviation of values, and Pooled SD = square-root of the sum of the squared within-participant standard deviations.

**Table 8. Summary of Results for Carbohydrates, %**

| Lab                     | A    | B     | C    | Mean | SD  | Method                                     |
|-------------------------|------|-------|------|------|-----|--------------------------------------------|
| E002                    | 7.02 | 6.95  | 7.16 | 7    | 0.1 | <i>Not specified</i>                       |
| E030                    | 18.9 | 19.8  | 20.2 | 19.6 | 0.7 | Calculation [100-(solids+protein+fat+ash)] |
| E033                    | 1.5  | 5.2   | 4.8  | 3.8  | 2   | Calculation [100-(solids+protein+fat+ash)] |
| E047                    | 11.7 | 12.21 |      | 12   | 0.4 | Calculation [100-(solids+protein+fat+ash)] |
| <b>N:</b>               |      |       |      | 4    |     |                                            |
| <b>Mean, Pooled SD:</b> |      |       |      | 10.5 | 2.2 |                                            |
| <b>SD:</b>              |      |       |      | 6.6  |     |                                            |

### 3.1.4. Ash

The ash values reported by laboratories participating in HAMQAP Exercise 5 are summarized in Table 9. The values from Lab E035 were implausibly large and therefore were omitted from the statistical analysis. The table also provides several summary values: *N* = number of values, Mean = mean of values, SD = standard deviation of values, and Pooled SD = square-root of the sum of the squared within-participant standard deviations.

**Table 9. Summary of Results for Ash, %**

| Lab                     | A    | B    | C    | Mean  | SD   | Method                            |
|-------------------------|------|------|------|-------|------|-----------------------------------|
| E002                    | 3.21 | 3.17 | 3.10 | 3.16  | 0.06 | Thermogravimetry                  |
| E009                    | 2.61 | 2.81 | 2.37 | 2.6   | 0.22 | <i>Not specified</i>              |
| E030                    | 2.69 | 2.67 | 2.67 | 2.68  | 0.01 | Muffle furnace                    |
| E033                    | 2.83 | 2.85 | 2.87 | 2.85  | 0.02 | Muffle furnace                    |
| E035                    | 84.6 | 78.4 | 69.9 | 77.63 | 7.38 | Muffle furnace                    |
| E047                    | 3.19 | 3.1  |      | 3.15  | 0.06 | Muffle furnace                    |
| <b>N:</b>               |      |      |      | 6     |      |                                   |
| <b>Mean, Pooled SD:</b> |      |      |      | 16.1  | 7.4  |                                   |
| <b>SD:</b>              |      |      |      | 29.5  |      |                                   |
| <b>N:</b>               |      |      |      | 5     |      | *removing outlying data from E035 |
| <b>Mean, Pooled SD:</b> |      |      |      | 2.87  | 0.24 |                                   |
| <b>SD:</b>              |      |      |      | 0.26  |      |                                   |

### 3.1.5. Solids

The solids values reported by laboratories participating in HAMQAP Exercise 5 are summarized in Table 10. The table also provides several summary values: *N* = number of values, Mean = mean of values, SD = standard deviation of values, and Pooled SD = square-root of the sum of the squared within-participant standard deviations.

**Table 10. Summary of Results for Solids, %**

| Lab                     | A     | B     | C     | Mean  | SD   | Method                    |
|-------------------------|-------|-------|-------|-------|------|---------------------------|
| E002                    | 97.67 | 97.7  | 97.73 | 97.7  | 0.03 | Thermogravimetry          |
| E030                    | 95.5  | 95.22 | 95.3  | 95.34 | 0.14 | Drying in Forced Air Oven |
| E033                    | 95.5  | 95.9  | 95.7  | 95.7  | 0.2  | Drying in Vacuum Oven     |
| E047                    | 96.3  | 96.5  |       | 96.4  | 0.14 | Drying in Forced Air Oven |
| <b>N:</b>               |       |       |       | 4     |      |                           |
| <b>Mean, Pooled SD:</b> |       |       |       | 96.3  | 0.3  |                           |
| <b>SD:</b>              |       |       |       | 1.0   |      |                           |

### 3.1.6. Calories

The calories values reported by laboratories participating in HAMQAP Exercise 5 are summarized in Table 11. The upper portion of the table summarizes the results as reported by the participants in the study. The lower portion of the table repeats this information, replacing the results for two laboratories who reported results in incorrect units. The table also provides several summary values: *N* = number of values, Mean = mean of values, SD = standard deviation of values, and Pooled SD = square-root of the sum of the squared within-participant standard deviations.

**Table 11. Summary of Results for Calories, kcal/100 g**

| HAMQAP Exercise 5 (as reported by participants) |        |        |        |       |       |                                                 |
|-------------------------------------------------|--------|--------|--------|-------|-------|-------------------------------------------------|
| Lab                                             | A      | B      | C      | Mean  | SD    | Method                                          |
| E002                                            | 655.44 | 653.26 | 658.93 | 656   | 3     | Calculation [9(fat)+4(protein)+4(carbohydrate)] |
| E030                                            | 0.629  | 0.62   | 0.619  | 0.623 | 0.006 | Calculation [9(fat)+4(protein)+4(carbohydrate)] |
| E033                                            | 682.4  | 670.2  | 667.3  | 673   | 8     | Calculation [9(fat)+4(protein)+4(carbohydrate)] |
| E047                                            | 0.642  | 0.639  |        | 0.641 | 0.002 | Calculation [9(fat)+4(protein)+4(carbohydrate)] |
| <b>N:</b>                                       |        |        |        | 4     |       |                                                 |
| <b>Mean, Pooled SD:</b>                         |        |        |        | 362.8 | 8.5   |                                                 |
| <b>SD:</b>                                      |        |        |        | 346.8 |       |                                                 |

| HAMQAP Exercise 5 (corrected for units) |        |        |        |       |      |                                                             |
|-----------------------------------------|--------|--------|--------|-------|------|-------------------------------------------------------------|
| Lab                                     | A      | B      | C      | Mean  | SD   | Method                                                      |
| E002                                    | 655.44 | 653.26 | 658.93 | 656   | 3    | Calculation [9(fat)+4(protein)+4(carbohydrate)]             |
| E030                                    | 629    | 620    | 619    | 623   | 6    | Calculation [9(fat)+4(protein)+4(carbohydrate)] - corrected |
| E033                                    | 682.4  | 670.2  | 667.3  | 673   | 8    | Calculation [9(fat)+4(protein)+4(carbohydrate)]             |
| E047                                    | 642    | 639    |        | 641   | 2    | Calculation [9(fat)+4(protein)+4(carbohydrate)] - corrected |
| <b>N:</b>                               |        |        |        | 4     |      |                                                             |
| <b>Mean, Pooled SD:</b>                 |        |        |        | 648.8 | 10.4 |                                                             |
| <b>SD:</b>                              |        |        |        | 20.9  |      |                                                             |

## 3.2. Trace Allergen Contaminants

### 3.2.1. Eurofins GeneScan

Data provided by Eurofins GeneScan for hazelnut protein in RM 8404 Almond Flour for Allergen Detection is reported in Table 12. No hazelnut protein was detected by the R-Biopharm ELISA kit above the lower limit of quantitation for all three samples.

**Table 12. Results provided by Eurofins for hazelnut allergen in RM 8404 using R-Biopharm Hazelnut Allergen ELISA**

| Sample Number | Result mg/kg (ppm) | LOQ mg/kg (ppm) |
|---------------|--------------------|-----------------|
| 1             | < 2.5              | 2.5             |
| 2             | < 2.5              | 2.5             |
| 3             | < 2.5              | 2.5             |

### 3.2.2. Food Allergen Research and Resource Program (FARRP) at the University of Nebraska Lincoln

Data provided by FARRP for hazelnut, peanut, pecan, walnut, and macadamia proteins in RM 8404 Almond Flour for Allergen Detection are reported in Table 13. No hazelnut, peanut, pecan, or walnut protein was detected by the Neogen and BioFront ELISA kits above the lower limits of quantitation for the assays in all three samples. A level of macadamia protein just above the LOQ was detected using the 3M kit. As described in Section **Error! Reference source not found.**, the likelihood of macadamia contamination of this almond flour is low based on low prevalence of macadamia nut production in the area where this material was sourced and is likely a result of potential cross-reactivity of the ELISA kit in the presence of the high level of almond protein.

**Table 13. Results provided by FARRP for various allergens in RM 8404**

| Allergen  | Test                                  | LOQ mg/kg (ppm) | Sample Result mg/kg (ppm) |       |       |
|-----------|---------------------------------------|-----------------|---------------------------|-------|-------|
|           |                                       |                 | 1                         | 2     | 3     |
| Hazelnut  | Neogen Veratox ELISA                  | 2.5             | < 2.5                     | < 2.5 | < 2.5 |
| Peanut    | Neogen Veratox ELISA                  | 2.5             | < 2.5                     | < 2.5 | < 2.5 |
| Pecan     | BioFront Technologies MonoTrace ELISA | 1.0             | < 1                       | < 1   | < 1   |
| Walnut    | BioFront Technologies MonoTrace ELISA | 1.0             | < 1                       | < 1   | < 1   |
| Macadamia | 3M ELISA                              | 0.33            | 0.37                      | 0.37  | 0.76  |

### 3.2.3. Microbac Laboratories

Data provided by Microbac Laboratories for hazelnut, peanut, pecan, walnut, pistachio, and cashew DNA in RM 8404 Almond Flour for Allergen Detection are reported in Table 14. No DNA from any tested source were detected above the method LOQ by the R-Biopharm PCR kits.

**Table 14. Results provided by Microbac Laboratories for various allergens in RM 8404 using R-Biopharm SureFood PCR Kits**

| Allergen  | LOQ<br>mg/kg<br>(ppm) | Qualitative Sample Result |          |          |
|-----------|-----------------------|---------------------------|----------|----------|
|           |                       | 1                         | 2        | 3        |
| Hazelnut  | 0.4                   | Negative                  | Negative | Negative |
| Peanut    | 0.4                   | Negative                  | Negative | Negative |
| Pecan     | 0.4                   | Negative                  | Negative | Negative |
| Walnut    | 0.4                   | Negative                  | Negative | Negative |
| Pistachio | 0.4                   | Negative                  | Negative | Negative |
| Cashew    | 0.4                   | Negative                  | Negative | Negative |

### 3.3. Customer Feedback

#### 3.3.1. Bia Diagnostics, LLC

As shown in Table 15, no cross-reactivity to or contamination with Brazil nut, cashew, coconut, hazelnut, macadamia, pecan, pistachio, or walnut was observed for RM 8404. The almond response for this kit was extremely high, with results above LOQ of 27 mg/kg (ppm) at a 1/62500 dilution, which equates to a result of over 1687500 mg/kg (ppm) almond protein. Because the almond antibody used in this kit was raised against roasted almonds, the blanched almonds used to prepare this RM may react more strongly than anticipated due to the abundance of intact, folded protein compared to those in roasted almonds as roasting may be more destructive to the protein 3D structure than blanching.

**Table 15. Summary of cross reactivity and contamination results provided by Bia Diagnostics ELISA testing of RM 8404.**

| Allergen   | Test LOQ<br>mg/kg (ppm) | Test Result<br>mg/kg (ppm) |
|------------|-------------------------|----------------------------|
| Almond     | 27                      | > 1 687 500                |
| Brazil Nut | 1                       | ND <sup>a</sup>            |
| Cashew     | 0.9                     | ND                         |
| Coconut    | 2                       | ND                         |
| Hazelnut   | 1                       | ND                         |
| Macadamia  | 0.3                     | ND                         |
| Pecan      | 0.66                    | ND                         |
| Pistachio  | 1                       | ND                         |
| Walnut     | 2                       | ND                         |

<sup>a</sup> ND = not detected

### 3.3.2. Eurofins Analytik GmbH

Eurofins Analytik GmbH reported that RM 8404 was homogenous and weighing and dissolution was straightforward to allow spiking experiments without issue. As shown in Table 16, no cross-reactivity to or contamination with bovine beta-lactoglobulin (a milk protein), casein (a milk protein), cashew, peanut, coconut, hazelnut, lupin, macadamia, Brazil nut, pecan, sesame, soy, egg, walnut, mustard, fish, oat, pistachio, celery, or wheat was observed for RM 8404 based on the ELISA and PCR tests conducted.

**Table 16. Results provided by Eurofins Analytik GmbH for RM 8404.**

| Allergen           | ELISA       |                | PCR             |                 |
|--------------------|-------------|----------------|-----------------|-----------------|
|                    | LOQ (mg/kg) | Result (mg/kg) | LOQ (mg/kg)     | Result (mg/kg)  |
| beta-lactoglobulin | 0.031       | < 0.031        | NT <sup>a</sup> | NT              |
| Brazil Nut         | 4           | < 4            | NT              | NT              |
| Casein             | 0.25        | < 0.25         | NT              | NT              |
| Cashew             | 2           | < 2            | Undefined       | ND <sup>b</sup> |
| Celery             | NT          | NT             | Undefined       | ND              |
| Coconut            | 2           | < 2            | NT              | NT              |
| Egg                | 0.31        | < 0.31         | NT              | NT              |
| Fish               | NT          | NT             | Undefined       | ND              |
| Hazelnut           | 0.16        | < 0.16         | Undefined       | ND              |
| Lupin              | 2           | < 2            | Undefined       | ND              |
| Macadamia Nut      | 1           | < 1            | NT              | NT              |
| Mustard            | 2           | < 2            | Undefined       | ND              |
| Oat                | NT          | NT             | Undefined       | ND              |
| Peanut             | 0.2         | < 0.2          | Undefined       | ND              |
| Pecan Nut          | 2           | < 2            | NT              | NT              |
| Pistachio          | NT          | NT             | Undefined       | ND              |
| Sesame             | 2           | < 2            | Undefined       | ND              |
| Soy                | 0.31        | < 0.31         | Undefined       | ND              |
| Walnut             | 2           | < 2            | Undefined       | ND              |
| Wheat              | NT          | NT             | Undefined       | ND              |

<sup>a</sup> NT = not tested

<sup>b</sup> ND = not detected

Additionally, Eurofins Analytik GmbH reported spiking RM 8404 into paprika and chocolate at 10 000 mg/kg (ppm) to evaluate recovery of their commodity assays. The paprika spiked with RM 8404 was tested with the almond ELISA and 8 400 mg/kg (ppm) of almond was recovered (84 %). The chocolate spiked with RM 8404 was tested with the almond ELISA and 7 700 mg/kg (ppm) of almond was recovered (77 %). Eurofins Analytik GmbH did not provide an acceptance range but indicated that recovery for all samples was in the range of normal acceptance criteria.

### 3.3.3. Eurofins Immunolab

Eurofins Immunolab tested two dilutions of RM 8404 and reported a mean activity of 95 % compared to the raw almond material applied for the calibration of the ELISA (Table 17).

**Table 17. Results for almond ELISA kit reactivity provided by Eurofins Immunolab for dilutions of RM 8404.**

| Kit-Lot | Sample   | Result mg/kg [ppm] |           | Reactivity [%] |
|---------|----------|--------------------|-----------|----------------|
|         |          | 1:100000           | 1:1000000 | Mean           |
| MDL-148 | Sample 1 | 9.60               | 0.80      | 88             |
|         | Sample 2 | 9.81               | 0.96      | 97             |
| MDL-149 | Sample 1 | 9.93               | 0.92      | 96             |
|         | Sample 2 | 9.75               | 0.98      | 98             |
| MDL-150 | Sample 1 | 9.87               | 0.99      | 99             |
|         | Sample 2 | 9.80               | 0.89      | 94             |
|         |          | Mean               |           | 95             |

### 3.3.4. Hygenia

Hygenia reported use of RM 8404 in quality control and cross reactivity studies and that the materials worked well, but that not enough material was provided for full characterization. No additional data or information was provided by Hygenia.

### 3.3.5. Neogen

Results for Neogen Veratox assay and Reveal 3D lateral flow kit testing of PBS and rice spiked with RM 8404 are provided in Table 18. Results for both assays are consistent with expectations.

**Table 18. Results provided by Neogen for ELISA and lateral flow kit testing of materials spiked with RM 8404.**

| Matrix     | Sample  | Veratox Almond ELISA |                    |          | Reveal 3D Almond Lateral Flow |          |
|------------|---------|----------------------|--------------------|----------|-------------------------------|----------|
|            |         | Spike mg/kg (ppm)    | Result mg/kg (ppm) | Recovery | Spike mg/kg (ppm)             | Result   |
| PBS        | Neg     | 0                    | 0.07               | -        | 0                             | Negative |
|            | RM 8404 | 5                    | 5.01               | 100 %    | 5                             | Positive |
|            |         | 10                   | 10.63              | 106 %    | 10                            | Positive |
|            |         | 20                   | 18.88              | 94 %     | 1000                          | Positive |
| Rice Flour | Neg     | 0                    | 0.18               | -        | 0                             | Negative |
|            | RM 8404 | 5                    | 5.32               | 106 %    | 5                             | Positive |
|            |         | 10                   | 10.42              | 104 %    | 10                            | Positive |
|            |         | 20                   | 23.16              | 116 %    | 1000                          | Positive |

## 4. Conclusions

### 4.1. Value Assignment for Proximates

As described in Section 2.3, available data for each measurand were used to provide an estimate of the mass fraction present in RM 8404 where  $x$  is the mean and  $U_{95}(x)$  is the 95 % confidence interval. The summary of these estimates is provided in Table 19, along with a summary of the datasets used to arrive at these estimates.

**Table 19. Summary of Estimates for Proximates in RM 8404**

| Analyte       | $x$    | $U_{95}(x)$ | Units      | Based on            |
|---------------|--------|-------------|------------|---------------------|
| Fat           | 56.02  | 8.03        | %          | HAMQAP              |
| Protein       | 23.29  | 1.18        | %          | HAMQAP, Eurofins    |
| Carbohydrates | 9.50   | 10.77       | %          | HAMQAP <sup>a</sup> |
| Ash           | 2.85   | 0.35        | %          | HAMQAP <sup>a</sup> |
| Solids        | 96.05  | 1.61        | %          | HAMQAP              |
| Calories      | 648.19 | 35.23       | kcal/100 g | HAMQAP              |

<sup>a</sup> Not all laboratories reported methods used.

## 4.2. Trace Allergen Contaminants

All testing results and customer feedback regarding the potential contamination of RM 8404 with other trace allergens is summarized in Table 20. Observed presence of macadamia nut is most likely explained as cross-reactivity of the 3M and Romer assays to the high levels of almond protein and therefore not of concern with respect to the proteins present in RM 8404.

**Table 20. Summary of all results for trace allergen contamination in RM 8404.**

| Allergen        | Presence<br>mg/kg<br>(ppm)                                     | Based on                                                                                                                                                                                              |
|-----------------|----------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| β-Lactoglobulin | < 0.031                                                        | Eurofins Analytik GmbH ELISA                                                                                                                                                                          |
| Brazil Nut      | < 1<br>< 1<br>< 4                                              | 3M ELISA<br>BioFront MonoTrace ELISA<br>Eurofins Analytik GmbH ELISA                                                                                                                                  |
| Casein          | < 0.25<br>< 0.4                                                | Eurofins Analytik GmbH ELISA<br>R-Biopharm SureFood PCR                                                                                                                                               |
| Cashew          | < 0.2<br>< 0.9<br>< 1<br>< 2<br>ND                             | Romer US-AgraQuant ELISA<br>3M ELISA<br>BioFront MonoTrace ELISA<br>Eurofins Analytik GmbH ELISA<br>Eurofins Analytik GmbH PCR                                                                        |
| Celery          | ND                                                             | Eurofins Analytik GmbH PCR                                                                                                                                                                            |
| Coconut         | < 2<br>< 2                                                     | 3M ELISA<br>Eurofins Analytik GmbH ELISA                                                                                                                                                              |
| Egg             | < 0.31                                                         | Eurofins Analytik GmbH ELISA                                                                                                                                                                          |
| Fish            | ND                                                             | Eurofins Analytik GmbH PCR                                                                                                                                                                            |
| Hazelnut        | < 0.16<br>< 0.3<br>< 0.4<br>< 1<br>< 1<br>< 2.5<br>< 2.5<br>ND | Eurofins Analytik GmbH ELISA<br>Romer US-AgraQuant ELISA<br>R-Biopharm SureFood PCR<br>BioFront MonoTrace ELISA<br>3M ELISA<br>R-Biopharm ELISA<br>Neogen Veratox ELISA<br>Eurofins Analytik GmbH PCR |
| Lupin           | < 2<br>ND                                                      | Eurofins Analytik GmbH ELISA<br>Eurofins Analytik GmbH PCR                                                                                                                                            |
| Macadamia       | < 0.3<br><b>0.5</b><br><b>35</b><br>< 1<br>< 2                 | 3M ELISA<br><b>3M ELISA</b><br><b>Romer US-AgraQuant ELISA</b><br>Eurofins Analytik GmbH ELISA<br>BioFront MonoTrace ELISA                                                                            |
| Mustard         | < 2<br>ND                                                      | Eurofins Analytik GmbH ELISA<br>Eurofins Analytik GmbH PCR                                                                                                                                            |

| Allergen  | Presence<br>mg/kg<br>(ppm)                    | Based on                                                                                                                                                              |
|-----------|-----------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Oat       | ND                                            | Eurofins Analytik GmbH PCR                                                                                                                                            |
| Peanut    | < 0.1<br>< 0.2<br>< 0.4<br>< 1<br>< 2.5<br>ND | Romer US-AgraQuant ELISA<br>Eurofins Analytik GmbH ELISA<br>R-Biopharm SureFood PCR<br>BioFront MonoTrace ELISA<br>Neogen Veratox ELISA<br>Eurofins Analytik GmbH PCR |
| Pecan     | < 0.4<br>< 0.66<br>< 1<br>< 1<br>< 2          | R-Biopharm SureFood PCR<br>3M ELISA<br>Romer US ELISA<br>BioFront MonoTrace ELISA<br>Eurofins Analytik GmbH ELISA                                                     |
| Pine Nut  | < 0.7<br>< 1                                  | Romer US ELISA<br>BioFront MonoTrace ELISA                                                                                                                            |
| Pistachio | < 0.1<br>< 0.4<br>< 1<br>< 1<br>ND            | Romer US-AgraQuant ELISA<br>R-Biopharm SureFood PCR<br>BioFront MonoTrace ELISA<br>3M ELISA<br>Eurofins Analytik GmbH PCR                                             |
| Sesame    | < 2<br>ND                                     | Eurofins Analytik GmbH ELISA<br>Eurofins Analytik GmbH PCR                                                                                                            |
| Soy       | < 0.31<br>ND                                  | Eurofins Analytik GmbH ELISA<br>Eurofins Analytik GmbH PCR                                                                                                            |
| Walnut    | < 0.3<br>< 0.4<br>< 1<br>< 2<br>< 2<br>ND     | Romer US-AgraQuant ELISA<br>R-Biopharm SureFood PCR<br>BioFront MonoTrace ELISA<br>3M ELISA<br>Eurofins Analytik GmbH ELISA<br>Eurofins Analytik GmbH PCR             |
| Wheat     | < 2<br>< 3<br>ND                              | Romer US-AgraQuant ELISA<br>Romer US-RIDASCREEN<br>Eurofins Analytik GmbH PCR                                                                                         |

<sup>a</sup> ND = not detected

## 4.3. Customer Feedback

Various laboratories spiked RM 8404 into various solutions or other food products to evaluate the recovery of their assays and fitness-for-purpose of the material. All assays responded acceptably and as expected.

## 5. References

- 1 Phillips MM, Bunk DM, Gonzalez CA, Jackson SA, Kauffman LR, Lippa KA, Rimmer CA, Semancik S, Winchester MR, Wood LJ (2020) Harnessing Measurement Science to Advance Food Safety. (National Institute of Standards and Technology, Gaithersburg, MD), NIST Special Publication (NIST SP) 1251. <https://doi.org/10.6028/NIST.SP.1251>
- 2 Barber CA, Burdette CQ, Hayes HV, Luvonga C, Phillips MM, Rimmer CA, Wood LJ, Yu LL (2020) NIST Health Assessment Measurements Quality Assurance Program: Exercise 5 Final Report. (National Institute of Standards and Technology, Gaithersburg, MD), NIST Internal Report (NIST IR) 8343. <https://nvlpubs.nist.gov/nistpubs/ir/2021/NIST.IR.8343.pdf>
- 3 AOAC INTERNATIONAL (2019) *Official Methods of Analysis 968.06 – Protein (Crude) in Animal Feed*. (AOAC INTERNATIONAL, Gaithersburg, MD).
- 4 Jones DB (1931) *Factors for Converting Percentages of Nitrogen in Food and Feeds into Percentages of Proteins* (United States Department of Agriculture, Beltsville, MD) Circular 183. [http://toolbox.foodcomp.info/References/Protein/Jones\\_1941%20nitrogen-protein%20conversion%20cir183.pdf](http://toolbox.foodcomp.info/References/Protein/Jones_1941%20nitrogen-protein%20conversion%20cir183.pdf) (accessed Oct 2021).
- 5 DerSimonian R, Laird N (1986) Meta-Analysis in Clinical Trials. *Controlled Clinical Trials* 7:177-188.
- 6 Joint Committee for Guides in Metrology (2008) *JCGM 100:2008 – Evaluation of Measurement Data — Guide to the Expression of Uncertainty in Measurement* (ISO GUM 1995 with Minor Corrections) <https://www.bipm.org/en/publications/guides>  
See also Taylor BN, Kuyatt CE (1994) Guidelines for Evaluating and Expressing the Uncertainty of NIST Measurement Results (National Institute of Standards and Technology, Gaithersburg, MD), NIST Technical Note (NIST TN) 1297. <https://doi.org/10.6028/NIST.TN.1297>
- 7 Joint Committee for Guides in Metrology (2008) *JCGM 101:2008 – Supplement 1 — Propagation of distributions using a Monte Carlo method*. <https://www.bipm.org/en/publications/guides>
- 8 Efron B, Tibshirani RJ (1993) *An Introduction to the Bootstrap* (Chapman & Hall, United Kingdom).
- 9 Searle S, Casella G, McCulloch C (1992) *Variance Components* (John Wiley, Hoboken, NJ).
- 10 Rukhin AL, Possolo A (2011) Laplace random effects models for interlaboratory studies. *Computational Statistics and Data Analysis* 55:1815-1827. <https://doi.org/10.1016/j.csda.2010.11.016>
- 11 Beauchamp CR, Camara JE, Carney J, Choquette SJ, Cole KD, DeRose PC, Duewer DL, Epstein MS, Kline MC, Lippa KA, Lucon E, Phinney KW, Polakoski M, Possolo A, Sharpless KE, Sieber JR, Toman B, Winchester MR, Windover D (2021) Metrological Tools for the Reference Materials and Reference Instruments of the NIST Materials Measurement Laboratory. (National Institute of Standards and Technology, Gaithersburg, MD), NIST Special Publication (NIST SP) 260-136, 2021 Edition. <https://doi.org/10.6028/NIST.SP.260-136-2021>

- 12 R-Biopharm AG (2021) RIDASCREEN FAST Hazelnut (R6802) Product Insert (R-Biopharm AG, Darmstadt, Germany) <https://food.r-biopharm.com/wp-content/uploads/2013/12/r6802-fast-hazelnut-21-03-12.pdf> (accessed Nov 2021).
- 13 Neogen Corporation (2018) Veratox for Hazelnut Allergen (8420) Product Insert (Neogen Corporation, Lansing, MI) <https://www.neogen.com/categories/allergens/veratox-hazelnut-allergen/> (accessed Nov 2021).
- 14 Neogen Corporation (2018) Veratox for Peanut Allergen (8430) Product Insert (Neogen Corporation, Lansing, MI) <https://www.neogen.com/categories/allergens/veratox-peanut-allergen/> (accessed Nov 2021).
- 15 BioFront Technologies MonoTrace ELISA Product Insert (BioFront Technologies, Tallahassee, FL) <https://www.biofronttech.com/category-brochure/3/> (accessed Nov 2021).
- 16 3M Food Safety (2017) Macadamia Protein ELISA Kit (E96MAC) Product Insert (3M Food Safety, St. Paul, MN) <https://multimedia.3m.com/mws/media/1622875O/3m-macadamia-protein-elisa-product-instructions.pdf> (accessed Nov 2021).
- 17 R-Biopharm AG (2020) SureFood Allergen Kits Product Insert (R-Biopharm AG, Darmstadt, Germany).  
[https://food.r-biopharm.com/wp-content/uploads/2020/08/surefood\\_allergen\\_kits\\_4p\\_en\\_2020-08.pdf](https://food.r-biopharm.com/wp-content/uploads/2020/08/surefood_allergen_kits_4p_en_2020-08.pdf) (accessed Nov 2021).

## Appendix A. Change Log

| Date   | Type of Edit | Change                                                                                                                                                                            | Location    |
|--------|--------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------|
| 9/2/26 | Editorial    | Updated formatting to meet current NIST template requirements                                                                                                                     | All         |
| 9/2/26 | Technical    | Additional detail added to describe the expansion of the uncertainty on the assigned value for protein based on the known reproducibility of the method used for value assignment | Section 3.3 |
| 9/2/26 | Technical    | Updated uncertainty value for protein                                                                                                                                             | Table 19    |