

Standard Reference Material[®] 2972b

25-Hydroxyvitamin D₂ and D₃ Calibration Solutions

CERTIFICATE OF ANALYSIS

Purpose: This Standard Reference Material (SRM) and its associated certified values are intended primarily for use in calibration of instruments and techniques used for the determination of the metabolites 25-hydroxyvitamin D₂ and 25-hydroxyvitamin D₃.

Description: A unit of SRM 2972b consists of fifteen two-milliliter ampoules, of three separate ethanolic solutions of individual vitamin D metabolites, five ampoules of each solution. Each two-milliliter ampoule contains approximately 1.2 mL of solution. The three separate solutions are described below:

25-Hydroxyvitamin D₃ Calibration Solution Level 1 [25(OH)D₃ Calibration Solution Level 1]
25-Hydroxyvitamin D₃ Calibration Solution Level 2 [25(OH)D₃ Calibration Solution Level 2]
25-Hydroxyvitamin D₂ Calibration Solution [25(OH)D₂ Calibration Solution]

Certified Values: The certified values for 25-hydroxyvitamin D₃ and 25-hydroxyvitamin D₂ in Table 1 are based on the analytical results determined using isotope dilution liquid chromatography with mass spectrometric detection (ID-LC-MS) and LC with absorbance detection (LC-absorbance). A NIST certified value is a value for which NIST has the highest confidence in its accuracy in that all known or suspected sources of bias have been investigated or taken into account [1]. The measurand is the total amount-of-substance of each analyte, expressed with either mass fraction (ng/g) or concentration units (nmol/L), in Table 1 [2]. Metrological traceability is to the International System of Units (SI) unit of mass expressed as mass fraction in nanograms analyte per gram of solution.

The certified amount-of-substance mass fraction values and concentration values for 25-hydroxyvitamin D₃ and 25-hydroxyvitamin D₂ are provided in Table 1. The certified amount-of-substance concentration values listed in Table 1 apply only to aliquots removed at 16 °C to 30 °C (see “Storage” and “Use”).

Additional Information: Additional information is provided in the appendices.

Period of Validity: The certified values delivered by **SRM 2972b** are valid within the measurement uncertainty specified until **31 August 2028**. The certified values are nullified if the material is stored or used improperly, damaged, contaminated, or otherwise modified.

Maintenance of Certified Values: NIST will monitor this SRM over the period of its validity. If substantive technical changes occur that affect the certification, NIST will issue an amended certificate through the NIST SRM website (<https://www.nist.gov/srm>) and notify registered users. SRM users can register online from a link available on the NIST SRM website or fill out the user registration form that is supplied with the SRM. Registration will facilitate notification. Before making use of any of the values delivered by this material, users should verify they have the most recent version of this documentation, available through the NIST SRM website (<https://www.nist.gov/srm>).

Carlos A. Gonzalez, Chief
Chemical Sciences Division

Steven J. Choquette, Director
Office of Reference Materials

Table 1. Certified Values for Vitamin D Metabolites in SRM 2972b

Vitamin D Metabolite	Mass Fraction ^(a) (ng/g)	Concentration ^(b) (nmol/L)
25-Hydroxyvitamin D ₃ (Level 1)	410.0 ± 14.9	806.2 ± 32.4
25-Hydroxyvitamin D ₃ (Level 2)	812.0 ± 29.2	1596.5 ± 64.1
25-Hydroxyvitamin D ₂	293.6 ± 9.1	560.4 ± 19.9

^(a) The uncertainty provided with each value is an expanded uncertainty about the weighted mean to cover the measurand with approximately 95 % confidence. The expanded uncertainty is calculated as $U = k u_c$, where u_c incorporates the observed difference between the results from the methods and their respective uncertainties, as well as uncertainties related to purity estimation and possible degradation of the solution over time, consistently with the ISO/JCGM Guide and with its Supplement 1, and k is the coverage factor ($k = 2$) corresponding to approximately 95 % confidence [3–7].

^(b) The amount-of-substance concentrations (nmol/L) were obtained by multiplying the certified values in mass fraction units by the density of ethanol (0.788 g/mL) at 22 °C and dividing by the relative molecular masses of 412.65 g/mol for 25-hydroxyvitamin D₂, and 400.64 g/mol for 25-hydroxyvitamin D₃. These concentrations are for use in the temperature range of 16 °C to 30 °C, and an allowance for the relative change in density over this temperature range of 0.87 % is incorporated in the uncertainties.

Safety: These solutions contain primarily ethanol, which is a flammable solvent. Open flames and sources of spark should be avoided while using this SRM. Use proper methods for disposal of flammable, potentially hazardous waste. Consult the Safety Data Sheet (SDS), enclosed with the SRM shipment, for health and safety information.

Storage: Sealed ampoules, as received, should be stored immediately in the dark at temperatures of –20 °C or below because of analyte instability at higher temperatures [8].

Use: Ampoules should be removed from the freezer and allowed to equilibrate to room temperature before weighing or volumetrically transferring. Due to the instability of the analytes in solution at temperatures greater than –20 °C, the total amount of time at room temperature for equilibrating and processing should be minimized to less than 3 h. Precautions should be taken to avoid exposure of ampoules and test portions to strong ultraviolet (UV) light and direct sunlight.

Test portions for use should be withdrawn immediately after opening the ampoules and should be processed or diluted without delay for the certified mass fraction values and concentration values to be valid within the stated uncertainty. Because of the volatility of ethanol, the certified values are NOT applicable to material stored in ampoules that have been opened for more than 2 min, even if they are resealed. The certified concentration values listed in Table 1 apply only to aliquots removed at 16 °C to 30 °C. If possible, samples should be placed in thermostatted compartments at 4 °C or colder during analysis.

Guidance for diluting the SRM 2972b solutions is provided in Appendix A.

Metrological Traceability

Metrological traceability of measurement results to a given reference must be established through an unbroken chain of calibrations and/or comparisons, each having stated uncertainties [9], using measurement standards that are appropriate for the property measured. The certified values in this calibration SRM are traceable to the International System of Units (SI) through such chains. These chains include: confirmation of the chemical identity and determination of the purity of primary standards, fitness evaluation of the solvent used to prepare each SRM solution, calibration of the devices used to determine mass, control of known influence factors such as temperature and UV radiation, and evaluation of the homogeneity and stability of each certified property as delivered in the SRM units.

Approaches to establishing the traceability of other measurement results of the certified property include calibration of a measurement process using the SRM as-is or by preparing in-house solutions through dilution of this SRM. The property values of in-house solutions can be made traceable to the SRM's certified value and through it to the SI by properly evaluating the uncertainties involved in procedures used in their preparation. Gravimetric and volumetric methods for preparing such in-house solutions are described in Appendix A Guidance for Diluting SRM 2972b 25-Hydroxyvitamin D₂ and D₃ Calibration Solutions. The property value uncertainties assigned to the in-house solutions must include the uncertainty of the SRM's certified value appropriately combined with the uncertainties of the preparative process.

The SI traceability of measurement results made using a measurement process calibrated to the SRM directly or to in-house solutions prepared from this SRM can then be established by properly evaluating the uncertainty of the calibration function resulting from the calibration process. The property value uncertainties assigned to measurement

results from the calibrated process must properly combine the calibration function uncertainty with the measurement process imprecision appropriate to the sample analyzed.

Guidance on evaluating and combining uncertainties is provided in reference 9.

REFERENCES

- [1] Beauchamp, C.R.; Camara, J.E.; Carney, J.; Choquette, S.J.; Cole, K.D.; DeRose, P.C.; Duewer, D.L.; Epstein, M.S.; Kline, M.C.; Lippa, K.A.; Lucon, E.; Molloy, J.; Nelson, M.A.; Phinney, K.W.; Polakoski, M.; Possolo, A.; Sander, L.C.; Schiel, J.E.; Sharpless, K.E.; Toman, B.; Winchester, M.R.; Windover, D.; *Metrological Tools for the Reference Materials and Reference Instruments of the NIST Material Measurement Laboratory*; NIST Special Publication (NIST SP) 260-136, 2021 edition; U.S. Government Printing Office: Washington, DC (2021); available at <https://nvlpubs.nist.gov/nistpubs/SpecialPublications/NIST.SP.260-136-2021.pdf> (accessed Oct 2022).
- [2] Thompson, A.; Taylor, B.N.; *Guide for the Use of the International System of Units (SI)*; NIST Special Publication 811; U.S. Government Printing Office: Washington, DC (2008); available at <https://www.nist.gov/pml/special-publication-811> (accessed Oct 2022).
- [3] JCGM 100:2008; *Evaluation of Measurement Data — Guide to the Expression of Uncertainty in Measurement*; (GUM 1995 with Minor Corrections), Joint Committee for Guides in Metrology (JCGM) (2008); available at <https://www.bipm.org/en/publications/guides> (accessed Oct 2022); see also Taylor, B.N.; Kuyatt, C.E.; *Guidelines for Evaluating and Expressing the Uncertainty of NIST Measurement Results*; NIST Technical Note 1297, U.S. Government Printing Office: Washington, DC (1994); available at <https://www.nist.gov/pml/nist-technical-note-1297> (accessed Oct 2022).
- [4] JCGM 101:2008; *Evaluation of Measurement Data — Supplement 1 to the “Guide to the Expression of Uncertainty in Measurement” — Propagation of Distributions using a Monte Carlo Method*; JCGM (2008); available at <https://www.bipm.org/en/publications/guides> (accessed Oct 2022).
- [5] DerSimonian, R.; Laird, N.; *Meta-Analysis in Clinical Trials*; Controlled Clin. Trials, Vol. 7, pp. 177–188 (1986).
- [6] Efron, B.; Tibshirani, R.J.; *An Introduction to the Bootstrap*; Chapman & Hall, London, UK (1993).
- [7] Searle, S.R.; Casella, G.; McCulloch, C.E.; *Variance Components*; John Wiley & Sons, New York, NY (1992).
- [8] Bedner, M.; Lippa, K.A.; *25-Hydroxyvitamin D Isomerizes to Pre-25-Hydroxyvitamin D in Solution: Considerations for Calibration in Clinical Measurements*; Analytical and Bioanalytical Chemistry, 407:8079-8086 (2015).
- [9] JCGM 200:2012; *International Vocabulary of Metrology - Basic and General Concepts and Associated Terms (VIM)* (2008 version with Minor Corrections), 3rd ed.; JCGM (2012); available at <https://www.bipm.org/en/publications/guides> (accessed Oct 2022).
- [10] EURACHEM; *EURACHEM/CITAC Guide: Quantifying Uncertainty in Analytical Measurement*; Ellison, S.; Williams, A. Eds.; 3rd ed. (2012); available at <https://www.eurachem.org/index.php/publications/guides> (accessed Oct 2022).

Certain commercial equipment, instruments, or materials may be identified in this Certificate of Analysis to adequately specify the experimental procedure. Such identification does not imply recommendation or endorsement by the National Institute of Standards and Technology, nor does it imply that the materials or equipment identified are necessarily the best available for the purpose.

Users of this SRM should ensure that the Certificate of Analysis in their possession is current. This can be accomplished by contacting the Office of Reference Materials 100 Bureau Drive, Stop 2300, Gaithersburg, MD 20899-2300; telephone (301) 975-2200; e-mail srminfo@nist.gov; or the Internet at <https://www.nist.gov/srm>.

***** End of Certificate of Analysis *****

APPENDIX A

Guidance for Diluting SRM 2972b 25-Hydroxyvitamin D₂ and D₃ Calibration Solutions

The ethanolic calibration solutions of SRM 2972b for 25-hydroxyvitamin D₃ and 25-hydroxyvitamin D₂ are of higher mass fraction and concentration than is typically encountered in human serum/plasma samples. Therefore, dilution of the calibration solutions may be required for analysis by many of the common vitamin D metabolite assays. Recommendations for dilution of the calibration solutions are as follows:

- 1) The solutions should be allowed to reach room temperature and be thoroughly mixed prior to opening the ampoules for dilution. However, due to the instability of the vitamin D metabolites at room temperature, care should be taken to minimize the total time at room temperature for equilibrating, diluting and analyzing to less than 3 h.
- 2) The most accurate results for dilution will be obtained using gravimetry with a calibrated analytical balance. Both the masses of the SRM 2972b solution and the total amount of solution after dilution are required to calculate the new mass fraction or concentration.
- 3) For assays that utilize volumetric measurements, use of either a gas-tight syringe or a positive displacement pipette (PDP) is recommended for solution transfer. If using a PDP, ensure all solution is delivered from the capillary (touching the tip of the capillary to the wall of the container may be required to fully deliver the correct volume). The best results for a volume dilution will be obtained if a volumetric flask is used to achieve the desired total volume.
- 4) If a positive displacement pipette is not available, an air-displacement pipette can be used. However, the errors in the amount of the ethanolic solution dispensed and the mass fraction or concentration of the diluted solution will be greater. Also, in both (3) and (4), use of pipettes that are out of calibration will result in additional error.
- 5) The choice of diluent does not matter, as long as the ethanol/analytes are soluble. Dilution with organic solvents such as alcohols or acetonitrile is preferable, but water can be used as the diluent to minimize solvent losses due to evaporation. All diluted SRM 2972b solutions should be stored in the dark in a sealed container (e.g., amber threaded bottle with a lined screw cap) to minimize mass fraction or concentration changes that could occur from evaporative losses. Solutions diluted with organic solvents should be stored at –20 °C until ready for analysis (up to 2 months) to minimize metabolite degradation. The viability of solutions that have been diluted with water, stored at –20 °C, and then equilibrated to room temperature for analysis has not been investigated at NIST. Therefore, it is recommended that samples diluted with water be analyzed without delay and discarded after 3 h.

* * * * * End of Appendix A * * * * *

APPENDIX B

Development of SRM 2972b was through collaboration between the National Institute of Standards and Technology (NIST) and the National Institutes of Health, Office of Dietary Supplements.

Coordination of the technical measurements leading to the certification of this SRM were performed by M. Bedner of the NIST Chemical Sciences Division. Analytical measurements at NIST were performed by M. Bedner, M.A. Nelson, and B.E. Lang of the NIST Chemical Sciences Division (CSD), as well as K.A. Lipka, formerly of CSD.

Statistical analysis was provided by J.H. Yen of the NIST Statistical Engineering Division.

Support aspects involved in the issuance of this SRM were coordinated through the NIST Office of Reference Materials.

Preparation and Analysis: The solutions were prepared gravimetrically at NIST from anhydrous ethanol and primary standards for 25-hydroxyvitamin D₂ obtained from IsoSciences (King of Prussia, PA) and 25-hydroxyvitamin D₃ obtained from the United States Pharmacopeia (Rockville, MD). The solutions of 25-hydroxyvitamin D₃ and 25-hydroxyvitamin D₂ were mixed overnight at ambient conditions (a minimum of 16 h). The morning following preparation, the solutions were chilled completely with ice and then aliquoted into 2 mL amber glass ampoules that had been purged with argon prior to addition of the solution. The ampoules were then flame-sealed. The masses of the primary standards and the total masses of the solutions were used to calculate the gravimetric mass fractions, which were used to confirm the analytical results described below.

NIST Analysis of Vitamin D Metabolites Using ID-LC-MS: Mass fractions of the 25-hydroxyvitamin D₂ and 25-hydroxyvitamin D₃ were measured at NIST using ID-LC-MS. Separate calibrants were prepared gravimetrically for each of the analytes at mass fractions intended to approximate the vitamin D metabolite levels in each solution of the SRM. Stable-isotope labeled internal standard solutions were used for each analyte and were added to the calibrants and the SRM 2972b samples. The internal standards consisted of separate solutions containing ²H₆-25-hydroxyvitamin D₃ and ²H₃-25-hydroxyvitamin D₂. Test portions ranged from approximately 200 mg to 500 mg, depending on the mass fraction of the metabolite being measured. For all SRM solutions, duplicate test portions from each of 10 ampoules selected using a stratified random sampling scheme were accurately weighed into 2-mL amber autosampler vials. An aliquot of the internal standard solution corresponding to the metabolite being measured was added, followed by mixing and injection. Details of the separation and a typical chromatogram for each of the two metabolites are provided in Figure B1. MS detection with selected ion monitoring was used for quantitation. 25-Hydroxyvitamin D₃ and ²H₆-25-hydroxyvitamin D₃ were monitored at *m/z* 383 and *m/z* 389, respectively. 25-Hydroxyvitamin D₂ and ²H₃-25-hydroxyvitamin D₂ were monitored at *m/z* 395 and *m/z* 398, respectively.

NIST Analysis of Vitamin D Metabolites Using LC-absorbance: Mass fraction values of the 25-hydroxyvitamin D₂ and 25-hydroxyvitamin D₃ were measured at NIST using LC-absorbance at 265 nm. Separate calibrants were prepared gravimetrically for each of the analytes at mass fractions intended to approximate the vitamin D metabolite levels in each solution of the SRM. The LC-absorbance measurements were calibrated using an external standard approach, and hence no internal standard was used. Aliquots from 10 ampoules, selected using a stratified random sampling scheme, were analyzed with LC-absorbance using both a C₁₈ and a pentafluorophenylpropyl (PFP) column. Representative chromatograms and the separation conditions for each of the metabolites are presented in Figure B2.

Purity Assessment: The mass fractions determined by gravimetric preparation, ID-LC-MS, and LC-absorbance were adjusted for the purity estimates of the primary standards, which were determined using LC-absorbance and three different stationary phases, LC-MS, thermogravimetric analysis, Karl Fischer titration, and quantitative nuclear magnetic resonance spectroscopy with an internal standard.

Homogeneity Assessment: The homogeneity of 25-hydroxyvitamin D₃ and 25-hydroxyvitamin D₂ in this SRM was assessed at NIST using the methods and test portion sizes described above, graphical analyses, and analyses of variance at the 5 % significance level. No significant inhomogeneity was observed for either of the vitamin D metabolites.

Value Assignment: The weighted mean of NIST results provided by ID-LC-MS and LC-absorbance were used to calculate assigned values for 25-hydroxyvitamin D₃ and 25-hydroxyvitamin D₂.

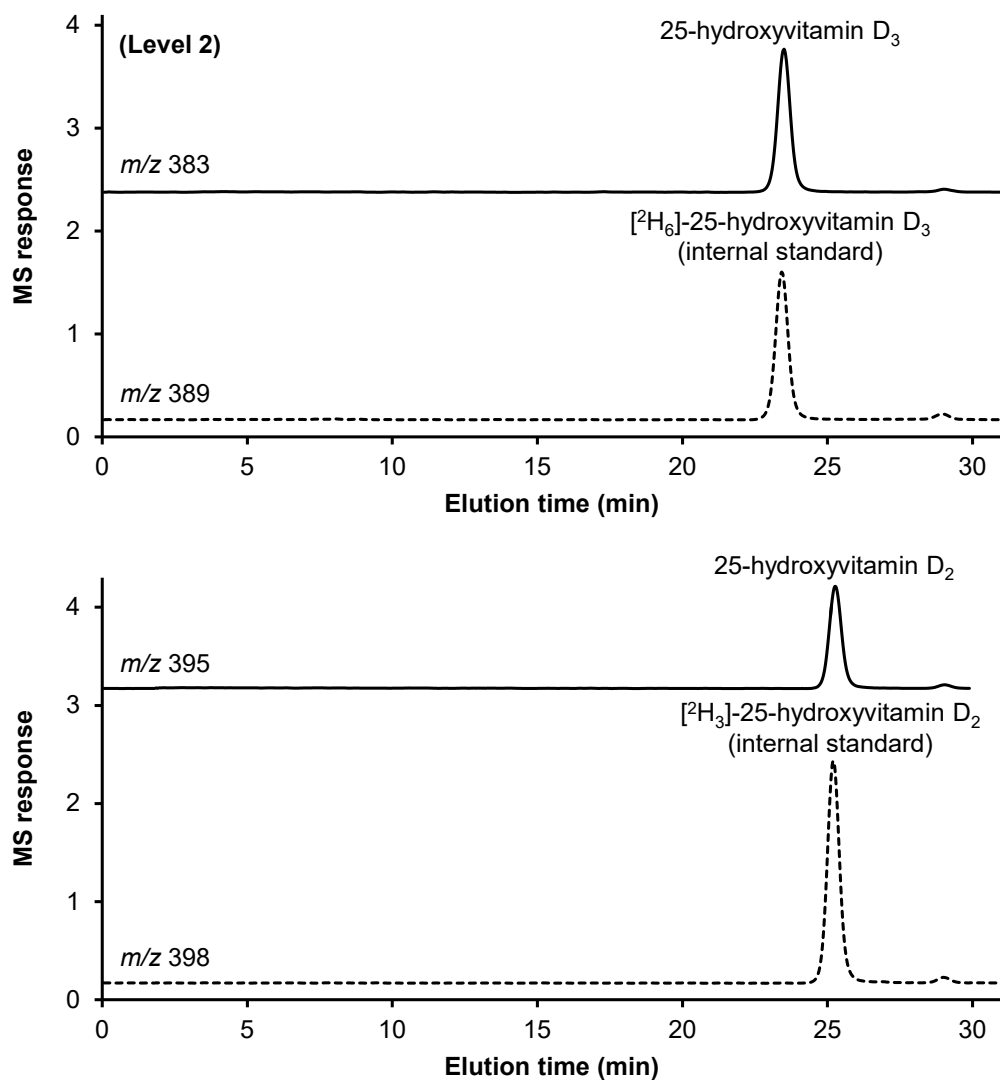


Figure B1. ID-LC-MS chromatograms of SRM 2972b 25-Hydroxyvitamin D₂ and D₃ Calibration Solutions. A pentafluorophenylpropyl column with dimensions of 150 mm × 4.6 mm ID and containing 2.7 µm diameter particles and an isocratic mobile phase of 78 % methanol, 22 % water were used with a column temperature of 15 °C. MS detection was achieved using selected ion monitoring at the indicated ions. The chromatograms for Level 1 of the 25-hydroxyvitamin D₃ solution are similar to those shown for Level 2.

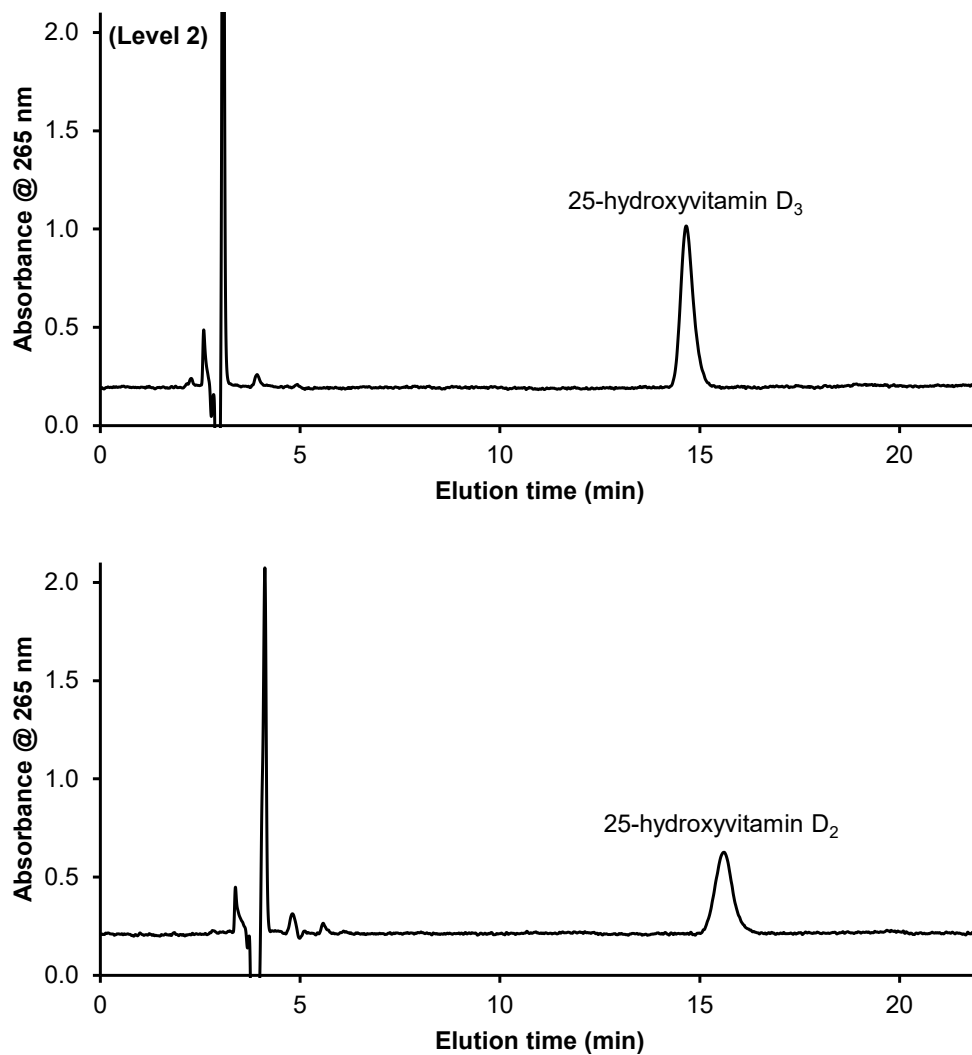


Figure B2. LC-Absorbance chromatograms of SRM 2972b 25-Hydroxyvitamin D₂ and D₃ Calibration Solutions. A C₁₈ column with dimensions of 250 mm × 4.6 mm ID and containing 5 μm particles and an isocratic mobile phase of 86 % methanol, 14 % water were used with a column temperature of 45 °C. Absorbance detection was achieved using 265 nm. The chromatograms for Level 1 of the 25-hydroxyvitamin D₃ solution are similar to those shown for Level 2.

***** End of Appendix B *****